

SMELL DISCETTES OLFACTION TEST

Instruction for Use

Statement

The instruction was originally written in German and is available in translation in other languages.

Note a more comprehensive version with video, in the customer section of the website. In addition to the brief instruction with three pictograms, on the back of the TestBox lid, detailed information is provided to interested users in this Instruction for use.



- Discettes must not be put in the mouth
- Keep the Discettes away from children
- Clarify whether the person being tested is allergic or hypersensitive to certain odours
- Dispose of unusable Discettes in accordance with the section below: «Return – Recycle – Dispose»

1 Intended Use

The 'Smell Discettes – Olfaction Test' (abbreviated: Smell-Test) is used for patient-individual screening to determine whether the patient's sense of smell is normal.

2 Patient groups and place of use

The Smell-Test is used in a clinical setting (hospital, medical practice, or pharmacy), on adults and children from the age of 6. The test requires the patient's cooperation.

3 Contra-Indications

The Discettes contain minimal amounts of odours that are commonly found in cosmetics and food products. Although we are not aware of any complications to date, patients should be asked before the test whether they have any allergies or intolerances to any of the eight odours. If the patient has an allergy to any of the odours used in the Smell-Test, the test must not be performed.

4 Smell-Test Applications

The Smell-Test is suitable for use in the following situations:

- Preoperative evaluation of the sense of smell prior to surgery on the nose and sinuses.
 - o Exclusion of a smell disorder that the patient is not aware of
 - o Preoperative documentation
 - o Postoperative documentation if necessary
- Presumptive Diagnosis of smell disorders.
 - o The Smell-Test can determine whether the sense of smell is normal
 - o If the Test shows that the sense of smell is not normal, further investigations must be carried out, for example, measurement of the olfactory threshold with N-butanol
- Clarification of any impairment of the sense of smell in diseases of the nose and paranasal sinuses.
 - o Acute rhinitis (flu, coronavirus infection)
 - o Chronic allergic rhinitis (hay fever, house dust mite allergy)

- o Acute and chronic sinusitis (nasal polyposis)
- Smell disorder after head/skull trauma
- Neurological diseases
- Smell disorder after contact with environmental toxins

5 Managing the Smell-Test Set

Before using the Smell-Test for the first time with patients, obtain a new complete SMELL Test-Set with the 8 different Discettes, in its TestBox. Our website helps you order the SMELL Test-Set in the “Product” section. The SMELL Test-Set will be supplied in an outer packaging, the shipping carton.

Art-#	GTIN	Name	Comments
SmellTest	7649988437519	SMELL Test-Set	The Smell-Test Set (treatment unit), consisting of the 8 different Discettes and documents, contained in the TestBox. Each Discette can be replaced individually.

6 Managing Smell-Test replacement Discettes

Depending on use, that is, how long each disc is left open, the smell may fade. Make sure each Discette is properly closed after use. If you notice the first signs of smell loss, obtain the corresponding replacement. It may be wise to purchase a full set of 8 replacement Discettes to reduce the risk of test failure: for example, in the event of sudden contamination. Our website helps you order replacement Discettes in the “Product” section.

On each Discette, you will find the manufacturing date for your reference: for example, 2025-09 (September of the year 2025).

Replacement Discettes for the SMELL Test-Set are supplied in a resealable (Gripp) bag and a shipping envelope.

Art-#	UDI-DI	Identification on Product	Remarks
Disc#01	7649988437526	1	Discette with the smell Coffee
Disc#02	7649988437533	2	Discette with the smell Vanilla
Disc#03	7649988437540	3	Discette with the smell Peach
Disc#04	7649988437557	4	Discette with the smell Grass
Disc#05	7649988437564	5	Discette with the smell Pineapple
Disc#06	7649988437571	6	Discette with the smell Rose
Disc#07	7649988437588	7	Discette with the smell Chocolate
Disc#08	7649988437595	8	Discette with the smell Fish

7 Smell-Test Execution

A Smell-Test takes about 15 minutes, including preparation and follow-up. When performed online, by a trained healthcare professional, it takes some minutes for a single patient. The general hygiene measures in force, in a clinical environment, must be observed. The sequential steps of the test are carried out online. The result is then available as a PDF document for further use. If preferred, a paper form is also available for copying in the test box. Further test forms can be downloaded, in the customer section of the website.

7.1 Smell-Test Preparation

To prepare the Smell Test for a patient, it is recommended to carefully wipe the TestBox and each Discette, without opening it, using a single-use disinfectant wipe.

- Wash and disinfect hands before testing to keep the materials clean.
- Defective Discettes, or Discettes without their typical odour, must no longer be used – replace them.
- The discs must always be handled carefully, without force, and must not be placed in the mouth.
- If a unilateral test is desired, the patient should close one nostril with a finger. In this case, the disc is presented to the patient by the healthcare professional.

7.2 Patient Preparation

To obtain a reliable result, no smoking, eating, or drinking should occur within 15 minutes before a test.

- Wash and disinfect hands before testing to keep the materials clean.
- People with severe allergies may not tolerate certain smells, which should be clarified with them beforehand.
- The discs must be handled carefully, without force, and must not be placed in the mouth.
- If a unilateral test is desired, one nostril can be closed with a finger.

7.3 Test a patient's sense of smell

Open the product's public website "www.smelldiscettes.com", by web-search or scanning the QR code on the back of the TestBox lid.


Select the language of your choice.

7.3.1 Access the OnlineTest

Click on SmellTest in the website menu and enter your login details to access the customer area of the website.



The login form is titled "Login". It contains three input fields: "Prüfer" (Name und vorname des Prüfer), "Klinik" (Klinik, in der der Test durchgeführt wird), and "Produktcode" (Code auf der Box vorhanden). A blue arrow points to a QR code icon next to the "Produktcode" field. Below the fields is a blue "SENDEN" button.

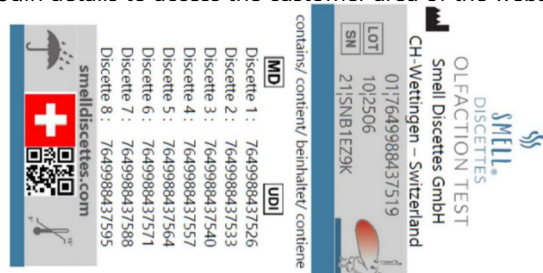

For this, click on the QR symbol in the last login line to scan the QR code (example above), next to the Swiss cross, on the label on the TestBox back or enter the 8-digit "21" serial number (SN), in case your device does not have a camera. After logging in, you are taken to the customer webpage, which contains information about usage and the OnlineTest. It is also possible to use copies of the paper form, for using the Smell-Test exceptionally without any digital tools.

7.3.2 Test Approach

1. The patient briefly opens a Discette to sniff at it without the Discette touching the nose.
2. The healthcare professional checks and wipes this opened Discette, presses it together with two fingers as shown in the picture on the right to securely close it, and returns it to the TestBox.
3. The patient must immediately choose the recognized smell from a selection of three illustrations (e.g., coffee – tomato – banana). One of the three illustrations closely matches the smell of the Discette. If the patient smells nothing, the most likely answer should be selected.
4. The cycle — open, smell, clean/close, and answer for each case — is repeated with all 8 Discettes. The order does not matter, but for clarity it is recommended to start with #1 and finish with Discette #8. If a disc has been contaminated by nasal secretions or other substances, it shall be replaced.
5. The OnlineTest immediately provides the result as a printable PDF document. The healthcare professional can also evaluate the personalized paper form immediately after the test's completion. Please refer to the "Evaluation" section below.
6. The physician discusses the result with the patient. To determine a possible extent of olfactory disorders, more in-depth testing methods are necessary e.g., determination of the olfactory threshold with n-Butanol.



7.4 Smell-Test Evaluation

The evaluation of a Smell-Test is usually carried out by the healthcare professional trained for this purpose, in the presence of the patient. Whether using the paper form or the online version, the Smell-Test result is calculated in the same way, as follows:

You get one point for each correct answer, so the result ranges from 0 to 8 points.

The correct answers are:

#1 coffee – #2 vanilla – #3 peach – #4 grass – #5 pineapple – #6 rose – #7 chocolate – #8 fish

Interpretation of the Smell-Test result

- 7 to 8 points: Normal olfactory ability
- 0 to 6 points: It is very likely that there is a disorder of the sense of smell (hyposmia, anosmia). It is also possible that the patient is not cooperating sufficiently.

7.5 Smell-Test Final Cleaning

To reduce the spread of germs and infections and thus contribute to maintaining an overall hygienic environment, the SMELL Test-Set should be integrated into the clinic's hygiene plan.

Each Discette and the TestBox must be wiped with a single-use disinfectant wipe after use.

The SMELL Test-Set with its Discettes is a non-critical product according to ISO 17664:2016-06, Annex C, and is not suitable for automated cleaning – like many stethoscopes.

7.6 Smell-Test Therapy Evaluation

If by the "SMELL Discettes Olfaction Test", the potential presence of an olfactory disorder was indicated, further examination may be necessary. In such an additional step, it is possible, for example, to determine the extent of the olfactory disorder or the olfactory threshold using a Test with n-butanol, which is not part of the Smell-Test.

8 Shelf-life, Storage and Replacement

The recommended shelf life of the SMELL Test-Set is generally 6 months, mainly for hygiene reasons. However, if the Smell-Test is used infrequently, well below 500 times, it may still provide reliable results beyond the six-month mark. Shortly before reaching 500 applications, the OnlineTest will indicate that replacement is needed. No later than one year after delivery, regardless of the number of uses, the Smell-Test must be fully replaced with a new SMELL Test-Set (TestBox with 8 different Discettes and documents). It is advisable to order a replacement early enough, ensuring uninterrupted testing.

A Smell-Test Set with tightly closed Discettes should be stored in a cool, dry place, ideally in a refrigerator at +2 to +15 °C. This helps slightly extend its lifespan, even with frequent patient testing.

9 Return – Recycle – Dispose

After its service life and final cleaning, a SMELL Test-Set should be returned to the distributor or directly to the manufacturer. This allows the discs to be separated for recycling of plastics. The same shipping packaging used for delivery can be used for this purpose. Alternatively, the discs may be sent to your local waste incineration facility, but not to a plastics recycling facility.

All paper and cardboard, on the other hand, may be recycled by yourself if you do not wish to return it.